

Unique X-Ray Visibility

## DERIVO<sup>®</sup> Embolisation Device



SMOOTH.  
RELIABLE.  
PRECISE.

**xcandis<sup>®</sup>**

ENGINEERING STROKE SOLUTIONS



## ADVANCED TECHNOLOGY FOR THE TREATMENT OF INTRACRANIAL ANEURYSMS

### SMOOTH.

- **Less friction during delivery and expansion** due to
  - BlueXide® Surface Finishing
  - Closed distal ends
- **Good vessel wall conformability** based on
  - Special braiding properties

### RELIABLE.

- **Secure wall apposition** because of
  - Flared ends
  - Closed distal ends
- **Better corrosion resistance and lower thrombogenicity** related to
  - BlueXide® Surface Finishing

### PRECISE.

- **Enhanced visualisation** as a result of
  - Nitinol Composite Wires with Platinum core
  - Three Platinum-Iridium X-Ray markers on both ends
- **Accurate placement** thanks to
  - Nitinol transport wire with s.e.c.u.r.e. GP Technology
  - Resheathability


## PROVEN TECHNOLOGY – SAFE AND EFFICIENT

### New composite wire concept for outstanding visibility of the DERIVO® contour

Treatment of left saccular ICA aneurysm with DERIVO® 5.0 mm x 20 mm



Excellent visibility of DERIVO® contour even in front of dense bone structures. View inside the lumen is possible.

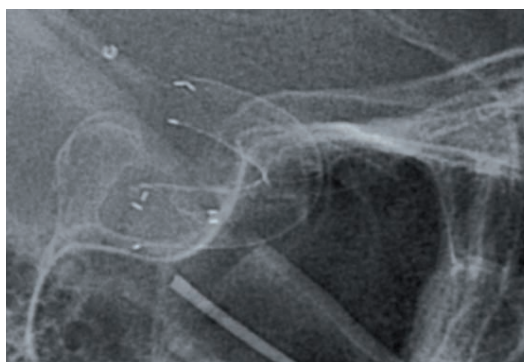


Opening of DERIVO® in tight curve is clearly visible.


*Images by courtesy of: Prof. Reith, Department of Neuroradiology, Saarland University Hospital, Homburg, Germany*

### Balanced mechanical properties for excellent clinical performance


Treatment of large right ICA aneurysm with DERIVO® 4.0 mm x 30 mm

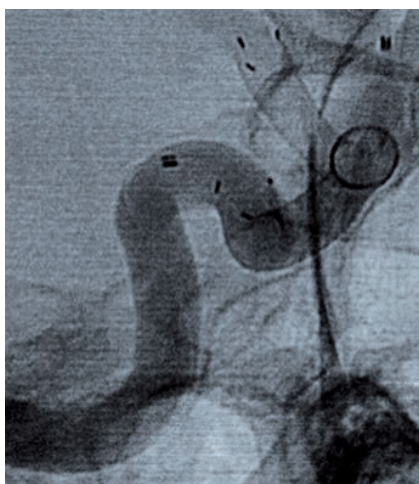

Perfect wall apposition: DERIVO® contour follows exactly the tortuous shape of the vessel.

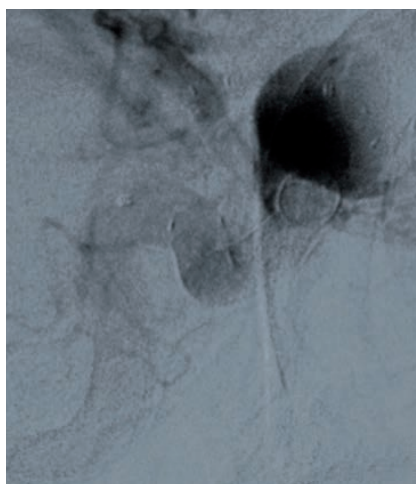


Immediate flow diversion effect after DERIVO® placement.

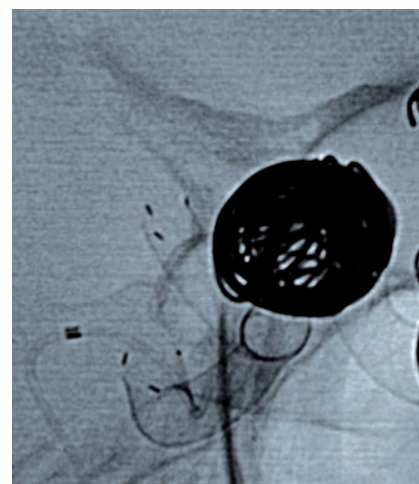

Excellent visibility of fully released DERIVO®.

*Images by courtesy of: Dr. Prothmann, Klinikum rechts der Isar, Department of Diagnostic and Interventional Neuroradiology, Technical University Munich, Germany*

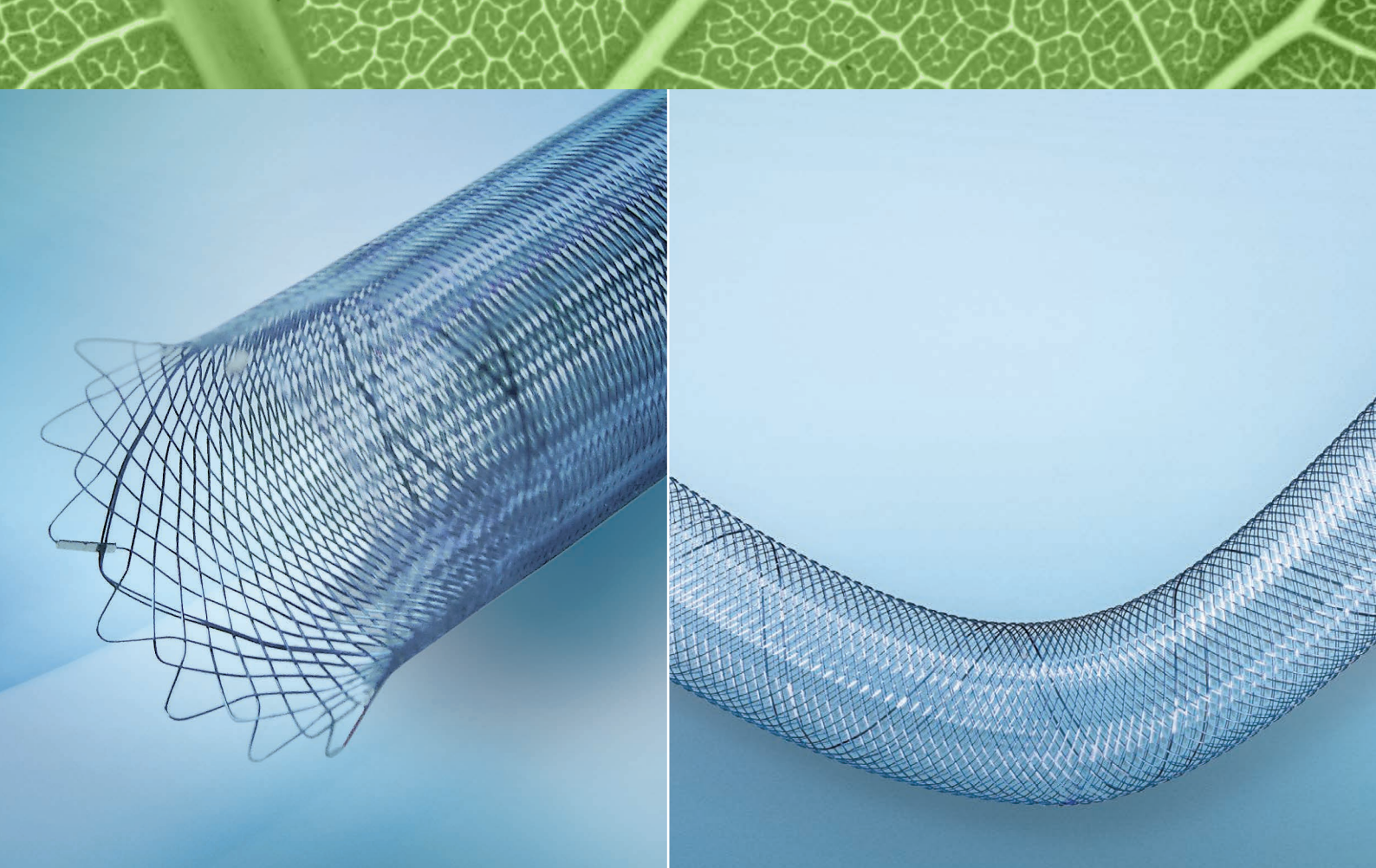

## FLOW – WHERE IT SHOULD BE

Acandis® is using the latest technological developments to ensure a smooth, reliable and precise treatment of intracranial aneurysms with the DERIVO® Embolisation Device.

### **BlueXide® Surface Finishing**

The Acandis® proprietary BlueXide® Surface Finishing Technology ensures less friction during delivery through the microcatheter as well as during expansion, making the opening of the device smooth and reliable. This finishing contributes to better corrosion resistance which might lead to lower thrombogenicity.

### **Nitinol Composite Wires**

The entire device consists of Nitinol Composite Wires with Platinum core leading to an outstanding visualisation of the contour and shape of the device under fluoroscopy.

### **X-Ray Markers**


Three Platinum-Iridium X-Ray markers are positioned on each end of the DERIVO® Embolisation Device for an accurate placement.

### **Closed Distal Ends**

The closed distal ends of the DERIVO® Embolisation Device help in delivering the device smoothly and releasing it simply, as they create less friction during the delivery through the microcatheter. Additionally these ends are less traumatic, even if the implant is oversized in the distal part of the vessel.

### **Flared Ends**

The DERIVO® Embolisation Device has flared ends for a secure wall apposition immediately after the initial distal opening, while the foreshortening on the proximal end is reduced.

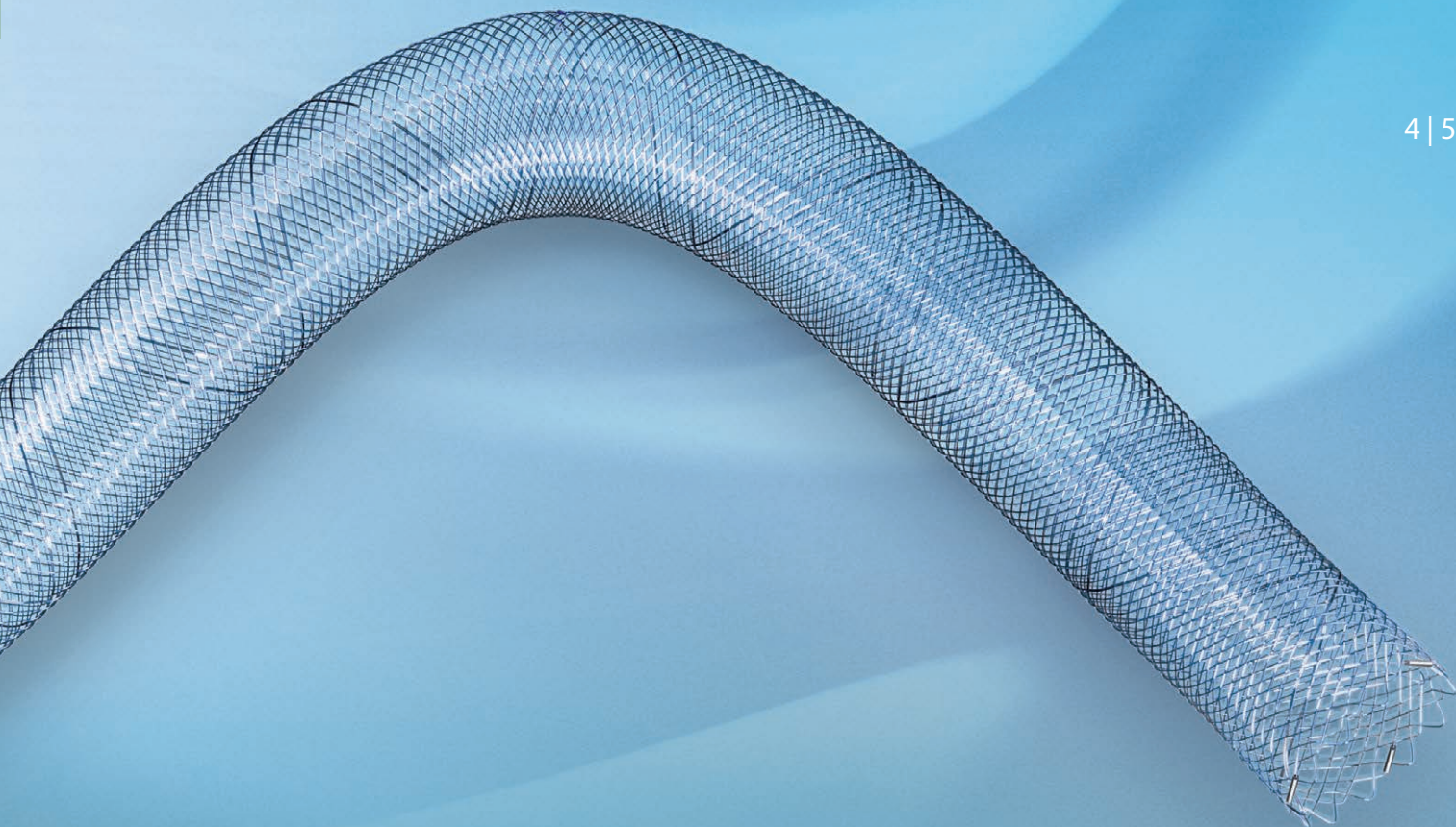

## Flow Diversion


The mesh density enables flow diversion away from the aneurysm while maintaining the flow into the side branches. Particle Image Velocimetry (PIV) proves the effectiveness of the DERIVO® Embolisation Device flow diversion properties.

## Vessel Wall Conformability



The braiding design ensures a good vessel wall conformability, even in highly variable vessel diameters and in tortuous anatomies.

## Velocity during Systole

Reference without Derivo® Embolisation Device

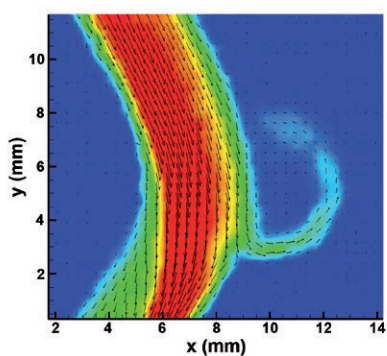


With Derivo® Embolisation Device

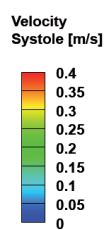
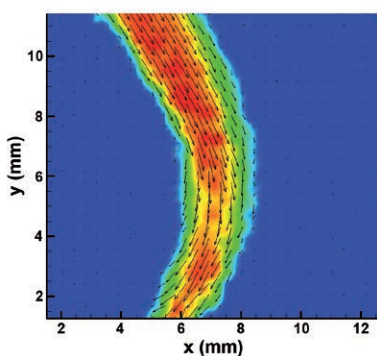


## Velocity during Diastole

Reference without Derivo® Embolisation Device

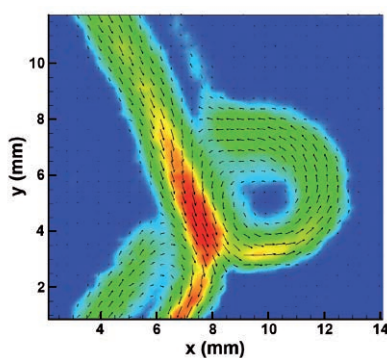

With Derivo® Embolisation Device

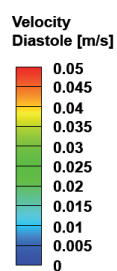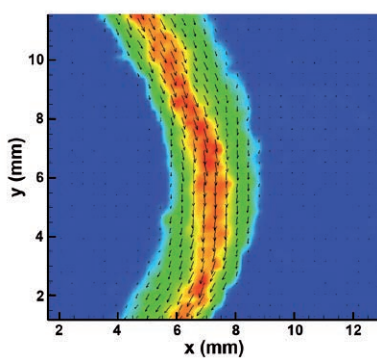

Particle Image Velocimetry (PIV) by courtesy of:  
Dept. of Cardiovascular Engineering RWTH  
Aachen (CVE/AME)

## PROCEDURE – RELIABLE AND EFFECTIVE

### s.e.c.u.r.e. GP Technology

The DERIVO® Embolisation Device is equipped with a Nitinol transport wire using the s.e.c.u.r.e. GP Technology engineered to meet the demands of a reliable and effective procedure.

**S-** safe

**E-** enhanced

**C-** controlled

**U-** unique

**R-** reliable


**E-** effective

The sleek surface of the transport wire changes into a unique – optically and tactile perceptible – checkered surface at the fluoroscopy marker point, to enhance the grip and push for a controlled and safe placement of the DERIVO® Embolisation Device.



### Resheathability

The device can be safely recaptured and repositioned if an adjustment and superior placement is needed.

### Tip Design

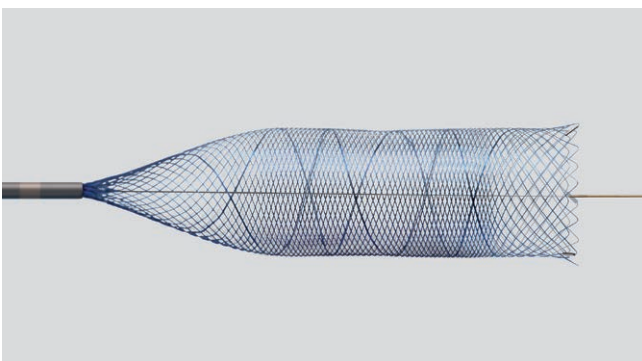

**With tip** – for additional distal support and retention of device access after release.

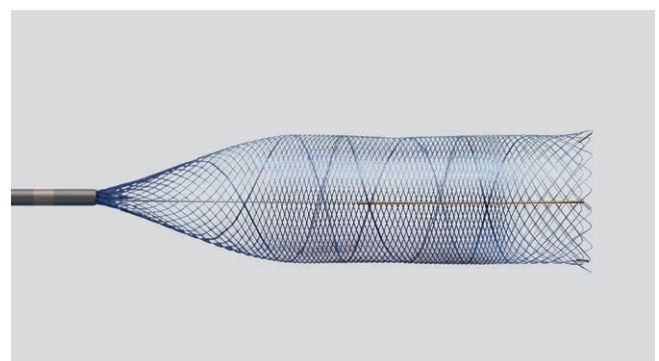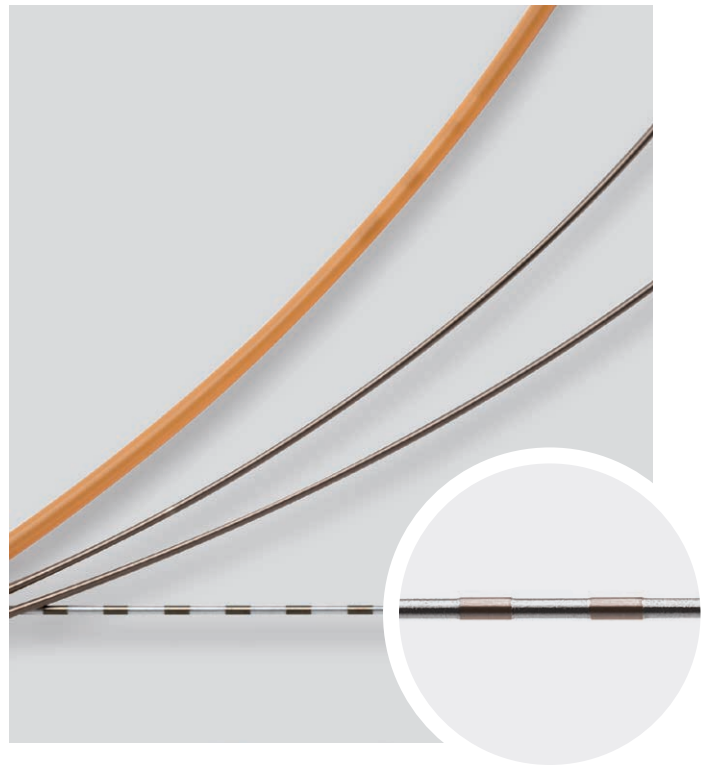

**Without tip (only applicable for 40 mm and 50 mm device lengths)** – for more flexibility and tip control in the treatment of long lesions.

## SIZING SUPPORT CHART – DERIVO® EMBOLISATION DEVICE

Labelled DERIVO® Dimensions (mm)	Reference Number		Unconstrained DERIVO® Dimensions (mm)	DERIVO® Lengths in corresponding Intended Use Diameters (mm)			
				Ø	3.5	3.0	2.5
		Ø	3.7	3.5	3.0	2.5	
3.5 × 15	01-000408	Device Length	10	15	20	25	
3.5 × 20	01-000409		13	20	27	32	
3.5 × 25	01-000410		16	25	35	41	
3.5 × 30	01-000411		19	30	41	48	
3.5 × 40	01-000415		25	40	53	66	
		Ø	4.2	4.0	3.5	3.0	
4.0 × 15	01-000381	Device Length	11	15	20	25	
4.0 × 20	01-000330		14	20	27	32	
4.0 × 25	01-000335		17	25	35	41	
4.0 × 30	01-000340		20	30	41	48	
4.0 × 40	01-000360		26	40	53	66	
		Ø	4.7	4.5	4.0	3.5	
4.5 × 15	01-000382	Device Length	11	15	20	25	
4.5 × 20	01-000331		14	20	27	32	
4.5 × 25	01-000336		17	25	35	41	
4.5 × 30	01-000341		20	30	41	48	
4.5 × 40	01-000361		26	40	53	66	
		Ø	5.2	5.0	4.5	4.0	
5.0 × 15	01-000383	Device Length	11	15	20	23	
5.0 × 20	01-000332		14	20	27	32	
5.0 × 25	01-000337		17	25	35	41	
5.0 × 30	01-000342		20	30	41	48	
5.0 × 40	01-000362		26	40	53	62	
5.0 × 50	01-000363	34	50	68	82		
		Ø	5.7	5.5	5.0	4.5	
5.5 × 15	01-000384	Device Length	11	15	20	23	
5.5 × 20	01-000333		14	20	27	32	
5.5 × 25	01-000338		17	25	35	41	
5.5 × 30	01-000343		20	30	41	48	
5.5 × 40	01-000364		26	40	53	62	
5.5 × 50	01-000365	34	50	68	82		
		Ø	6.2	6.0	5.5	5.0	
6.0 × 15	01-000385	Device Length	11	15	20	23	
6.0 × 20	01-000334		14	20	27	32	
6.0 × 25	01-000339		17	25	35	41	
6.0 × 30	01-000344		20	30	41	48	
6.0 × 40	01-000366		26	40	53	62	
6.0 × 50	01-000367	34	50	68	82		

Note: all indicated lengths can vary within a tolerance range of +/- 1mm

# ORDERING INFORMATION

Labelled DERIVO® Diameter (mm)	Labelled DERIVO® Length (mm)	Reference Number	Recommended Vessel Diameter (mm)	Required Microcatheter for Delivery ** (inch)
3.5	15	01-000408	2.5 – 3.5	0.027
	20	01-000409		
	25	01-000410		
	30	01-000411		
	40	01-000415*		
4.0	15	01-000381	3.0 – 4.0	
	20	01-000330		
	25	01-000335		
	30	01-000340		
	40	01-000360*		
4.5	15	01-000382	3.5 – 4.5	
	20	01-000331		
	25	01-000336		
	30	01-000341		
	40	01-000361*		
5.0	15	01-000383	4.0 – 5.0	
	20	01-000332		
	25	01-000337		
	30	01-000342		
	40	01-000362*		
	50	01-000363*		
5.5	15	01-000384	4.5 – 5.5	
	20	01-000333		
	25	01-000338		
	30	01-000343		
	40	01-000364*		
	50	01-000365*		
6.0	15	01-000385	5.0 – 6.0	
	20	01-000334		
	25	01-000339		
	30	01-000344		
	40	01-000366*		
	50	01-000367*		

\* Indicated on package as „without Tip“ as the tip always stays inside the stent for the 40 mm and 50 mm length

\*\* Please contact your local Acandis® representative for information on compatible microcatheters

All changes or modifications, may they be technical or other, or changes in the availability of products are expressly reserved.

Distributed by:



**ACANDIS GmbH & Co. KG**  
Theodor-Fahrner-Str. 6  
75177 Pforzheim  
Germany



Tel: +49 7231 155 00 0  
Fax: +49 7231 155 00 129  
E-Mail: info@acandis.com  
[www.acandis.com](http://www.acandis.com)